

A dentist in light blue scrubs and a white surgical mask is performing a procedure on a patient's teeth. The patient is lying back in a dental chair. The dentist is wearing blue gloves and using a dental instrument. The background is a light green wall with a large, stylized graphic of a dental arch in blue and white.

Aragó Sutures

PTFE

Polytetrafluoroethylene






A premium suture
for maxillofacial surgery
professionals






Laboratorio Aragón

Salvador Mundí, 11
08017 Barcelona
Tel. 93 204 22 12
Fax 93 203 26 09
export@laboratorioarago.com
www.laboratorioarago.com


3/8 circle needle cutting

| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|------|--------|---|----------------------|----------|-----|--------|--------------|
| TB8 | 11 mm |  | 75 cm | Incoloro | 6/0 | 0,7 | 56001 |
| | | | 75 cm | | 5/0 | 1 | 56002 |
| TB10 | 13 mm |  | 75 cm | Incoloro | 6/0 | 0,7 | 56003 |
| | | | 75 cm | | 5/0 | 1 | 56004 |
| | | | 75 cm | | 4/0 | 1,5 | 56005 |
| TB12 | 16 mm |  | 75 cm | Incoloro | 6/0 | 0,7 | 56006 |
| | | | 75 cm | | 5/0 | 1 | 56007 |
| | | | 75 cm | | 4/0 | 1,5 | 56008 |
| | | | 75 cm | | 3/0 | 2 | 56009 |
| TB15 | 19 mm |  | 75 cm | Incoloro | 5/0 | 1 | 56010 |
| | | | 75 cm | | 4/0 | 1,5 | 56011 |
| | | | 75 cm | | 3/0 | 2 | 56012 |
| | | | 75 cm | | 2/0 | 3 | 56013 |
| TB19 | 24 mm |  | 75 cm | Incoloro | 3/0 | 2 | 56014 |
| | | | 75 cm | | 2/0 | 3 | 56015 |

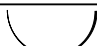
1/2 circle needle cutting

| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|------|--------|---|----------------------|----------|-----|--------|--------------|
| TC8 | 12 mm |  | 75 cm | Incoloro | 6/0 | 0,7 | 56016 |
| | | | 75 cm | | 5/0 | 1 | 56017 |
| TC12 | 20 mm |  | 75 cm | Incoloro | 4/0 | 1,5 | 56024 |
| | | | 75 cm | | 3/0 | 2 | 56025 |
| TC14 | 22 mm |  | 75 cm | Incoloro | 4/0 | 1,5 | 56026 |
| | | | 75 cm | | 3/0 | 2 | 56027 |
| | | | 75 cm | | 2/0 | 3 | 56028 |

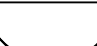


3/8 circle round needle trocar point

| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|-------|--------|---|----------------------|----------|-----|--------|--------------|
| PTB13 | 18 mm |  | 75 cm | Incoloro | 4/0 | 1,5 | 56018 |
| | | | 75 cm | | 3/0 | 2 | 56019 |



1/2 circle round needle trocar point


| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|------|--------|---|----------------------|----------|-----|--------|--------------|
| PTC8 | 10 mm |  | 75 cm | Incoloro | 4/0 | 1,5 | 56020 |

3/8 circle cutting needle precision point

| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|-------|--------|---|----------------------|----------|-----|--------|--------------|
| ATB8 | 11 mm |  | 45 cm | Incoloro | 7/0 | 0,5 | 56021 |
| | | | 45 cm | | 6/0 | 0,7 | 56022 |
| ATB10 | 13 mm |  | 45 cm | Incoloro | 5/0 | 1 | 56023 |
| ATB12 | 15 mm |  | 45 cm | Incoloro | 4/0 | 1,5 | 56029 |
| | | | 45 cm | | 3/0 | 2 | 56030 |

3/8 circle cutting needle black cutting point

| REF | NEEDLE | | LENGTH OF THE THREAD | COLOUR | USP | METRIC | CODE |
|-------|--------|---|----------------------|----------|-----|--------|--------------|
| BTB10 | 13 mm |  | 75 cm | Incoloro | 5/0 | 1 | 56031 |
| | | | 75 cm | | 4/0 | 1,5 | 56032 |
| BTB12 | 16 mm |  | 75 cm | Incoloro | 4/0 | 1,5 | 56033 |
| | | | 75 cm | | 3/0 | 2 | 56034 |
| | | | 75 cm | | 2/0 | 3 | 56035 |










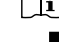
 Packaging in boxes of 12 units


 Check available variations in thread length and needle with Laboratorio Aragón

PTFE

Extreme Softness

| | |
|------------------------|--|
| Structure | Monofilament |
| Composition | Polytetrafluoroethylene (PTFE) |
| Colour | Incoloro |
| Sterilization | Ethylene oxide |
| Absorption | Non absorbable |
| Characteristics | <ul style="list-style-type: none"> • High chemical inertness and no tissue reaction • High impermeability • Extreme softness, it has a very low coefficient of friction • Maintains its qualities in moist tissues • Minimizes damage around the suture space preventing bacterial colonization in the deeper layers of the wound • Compared to classic monofilament sutures, it provides a durable and safe knotting of the thread. The ends of the suture do not cause irritation of the cheek, lips and tongue unlike materials based on other monofilament fibres that can additionally contribute to the formation of changes in the mucosa • In contrast to sutures with braided filaments, especially those that are made of absorbable materials, PTFE does not absorb blood, saliva, bacteria or food debris • Promotes wound healing |
| Indications | <ul style="list-style-type: none"> • Implantology • Bone reconstruction • Periodontology • Minimally invasive bone augmentation: the use of PTFE provides a significant reduction in postoperative inflammatory complications as well as the best and fastest wound healing • In open-healing approach • In retention suture applied to wounds on the skin and mucous membrane |
| Packaging | Boxes of 12 units |

-  Nº lote
-  Referencia catálogo
-  Un solo uso
-  No reesterilizar
-  Caducidad yyyy-mm año-mes
-  Esterilizado por óxido de etileno
-  Preservar de la luz de sol
-  No utilizar si el envase unitario está abierto o deteriorado
-  Ver instrucciones
-  Fabricante

 Marcado CE y nº de identificación del Organismo Notificado. El producto cumple con los Requisitos Esenciales de la Directiva de Productos Sanitarios 93/42/CEE.

The Premium suture

Polytetrafluoroethylene (PTFE)



For maxillofacial surgery professionals